


tobii dynavox TD Phone Connection Instructions

[Home](#) » [tobii dynavox](#) » tobii dynavox TD Phone Connection Instructions 

tobii dynavox

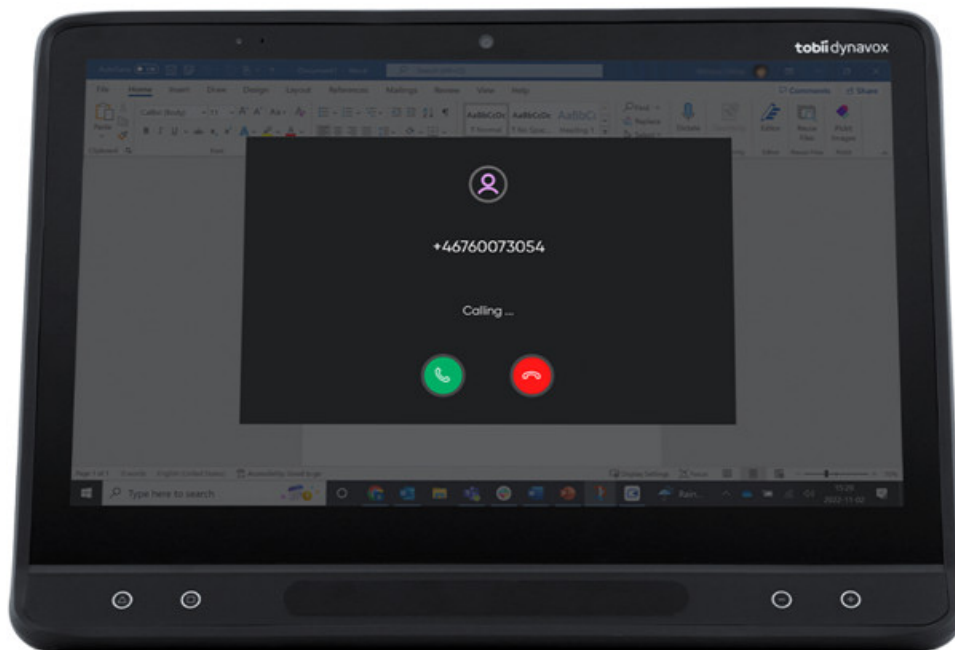

**TD Phone
Connection Guide**

Contents

- 1 TD Phone Connection
- 2 Before Connecting
- 3 Connect your Mobile Phone
- 4 Tour of TD Phone
- 5 Using TD Phone
- 6 Send a Text Message
- 7 6 page start
 - 7.1 Documents / Resources
 - 7.1.1 References

TD Phone Connection



TD Phone app

TD Phone allows you to text and interact with your smart phone through eye gaze on your Tobii Dynavox I-Series device.

Before Connecting

Make sure that TD Phone is installed on your Tobii Dynavox I-13 or I-16 device. To download the TD Phone installer, visit www.tobiidynavox.com/pages/td-phone-product-support?tab=1 and select TD Phone Software Bundle Installer.



Note

If you have a closed I-13 or I-16 device, use another computer to download the TD Phone installer from the link above. Use a USB storage device to move the TD Phone installer file from the computer to your Tobii Dynavox device.

Initial setup requires input using the phone touch screen. If you are not able to access the phone touch screen, ask an assistant to complete these steps for you.

Connect your Mobile Phone

The first time that you launch TD Phone, the on-screen connection assistant will start automatically. Follow the instructions carefully in each step to connect your mobile phone and start using TD Phone. You can open the connection assistant again at any time:


1. Look below the screen to open the Off-Screen Menu.
2. Select Dashboard.
3. Select Settings.
4. Select the right arrow.
5. Select My Phone.
6. Select the Connect button in the Connect Mobile Phone section.



Essential for Successful Connection

Your mobile phone must always remain unlocked throughout the connection process.

Tour of TD Phone

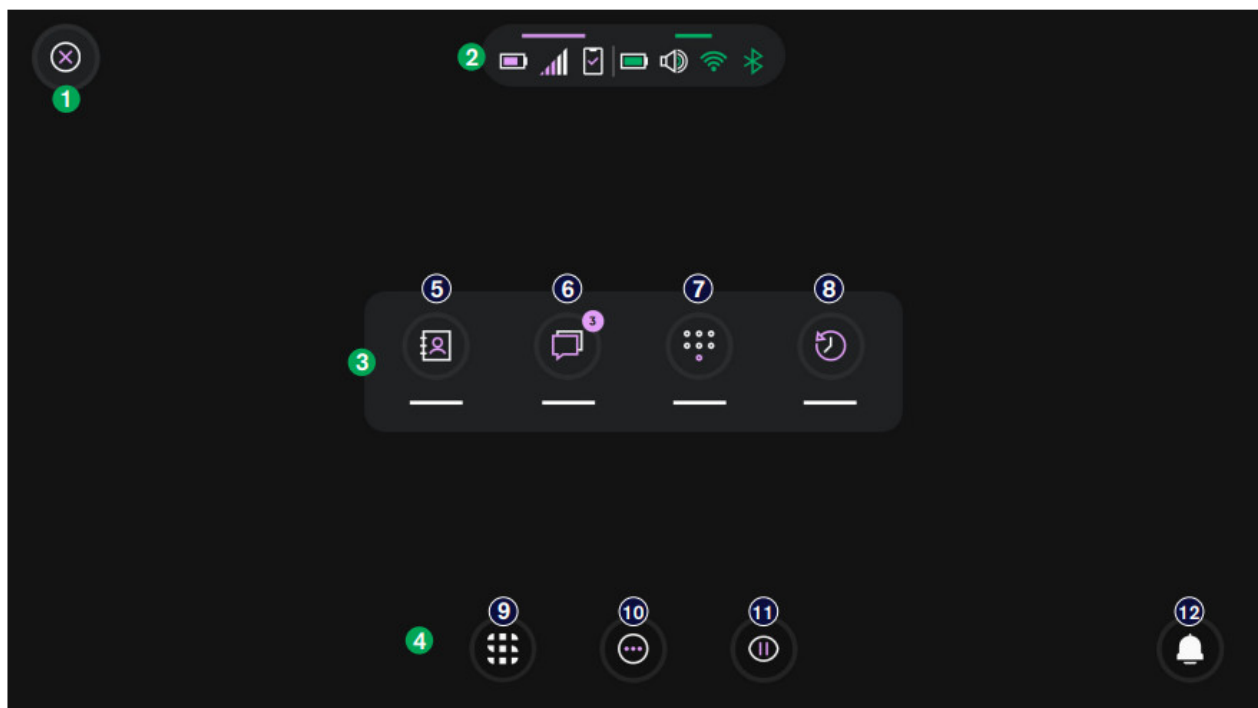

Main Page

1. Open/close main menu
2. Status Bar
3. Main Menu
4. Off-Screen Menu:
Appears when you look below the eye tracker

Main Menu

5. Contacts:
Browse your contacts. Select a contact to message or call them.
6. Messages:
View, send, and receive text messages

7. Keypad:

Call or message any phone number

8. Recents:

View your recent incoming and outgoing calls

Off-Screen Menu

9. Switcher:

Move between apps.

10. Dashboard:


TD Phone Settings.

11. Pause/Resume eye tracking:

Stop and start eye tracking

12. Do Not Disturb:

Turn off TD Phone notifications

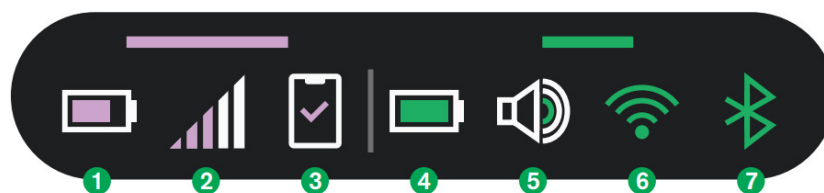

Status Bar indicators

1. Mobile phone battery

2. Mobile phone signal


3. Mobile phone connection to device

4. Device battery

5. Device volume level

6. Device Wifi signal

7. Device Bluetooth status

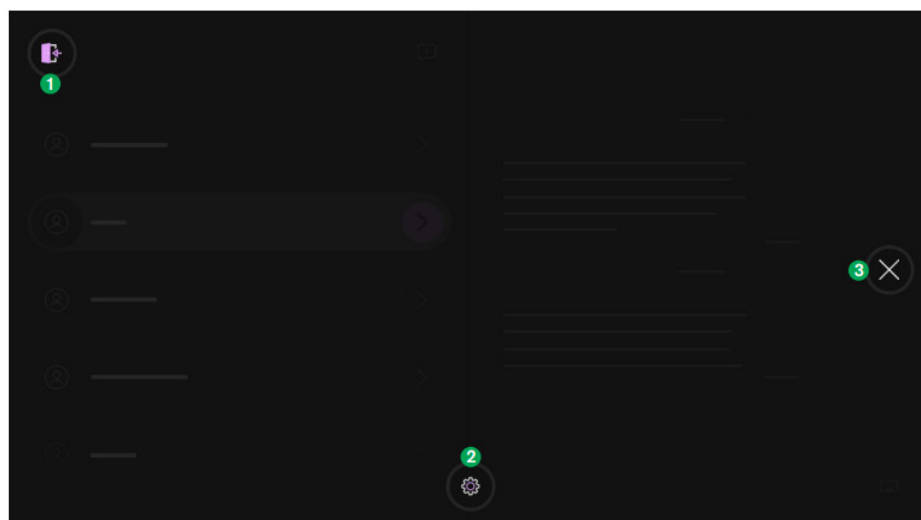

Dashboard

1. Exit TD Phone

2. Settings:

Where you connect your mobile phone to TD Phone and access other settings such as dwell time.

3. Close Dashboard:

Returns to the Main Menu

Using TD Phone




Note







- TD Phone requires that your mobile phone remain on. It does not need to remain unlocked.
- Your mobile phone must remain within 3-5 meters (10-16 feet) of your I-Series device to maintain a good Bluetooth connection. The actual range varies, depending on your mobile phone model.







Send a Text Message




—

6 page start

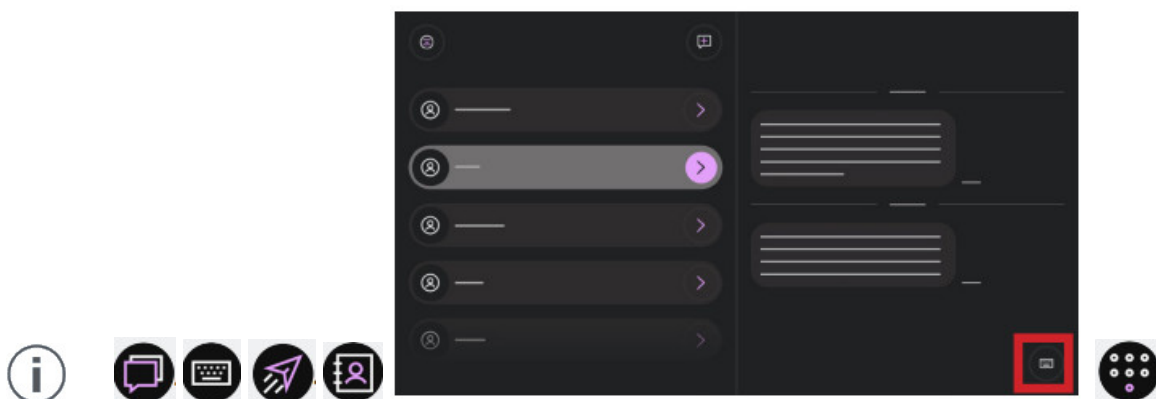

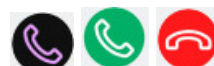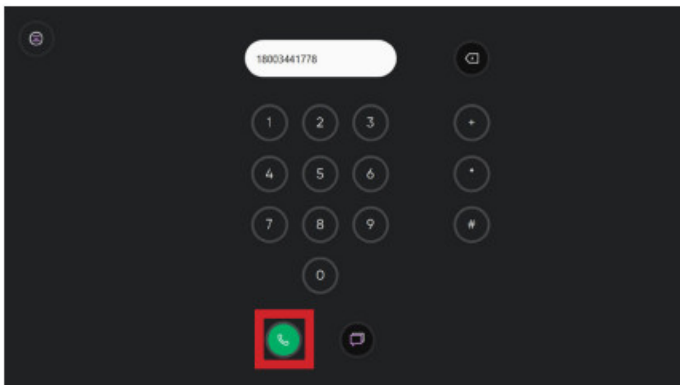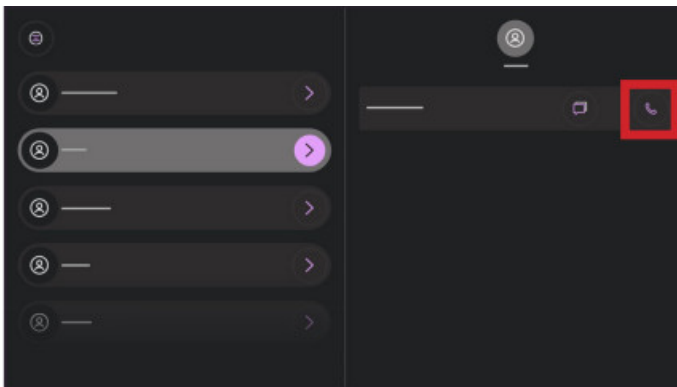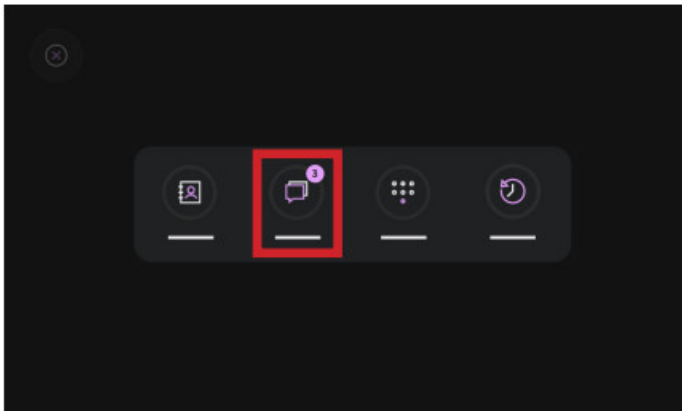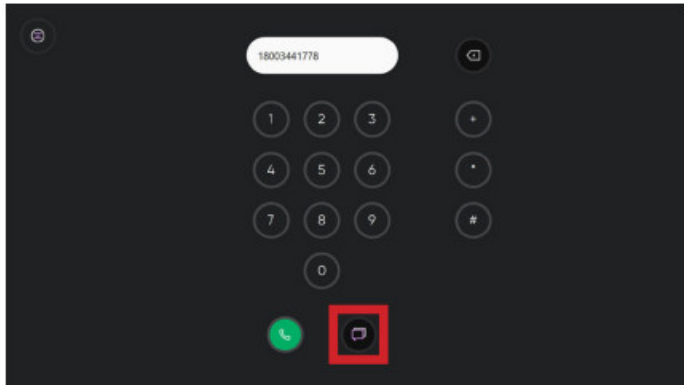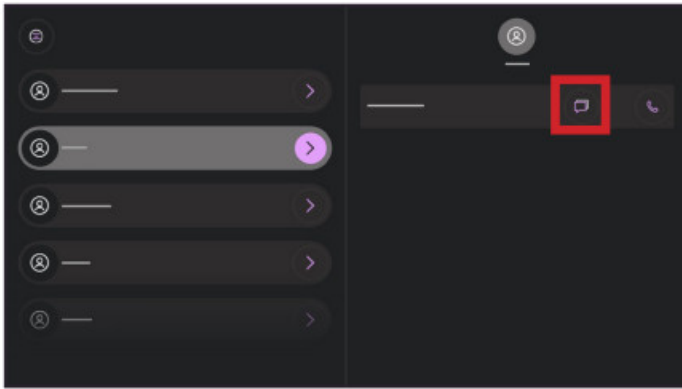

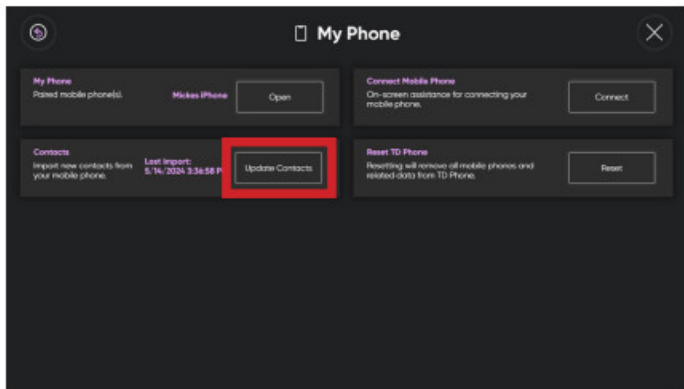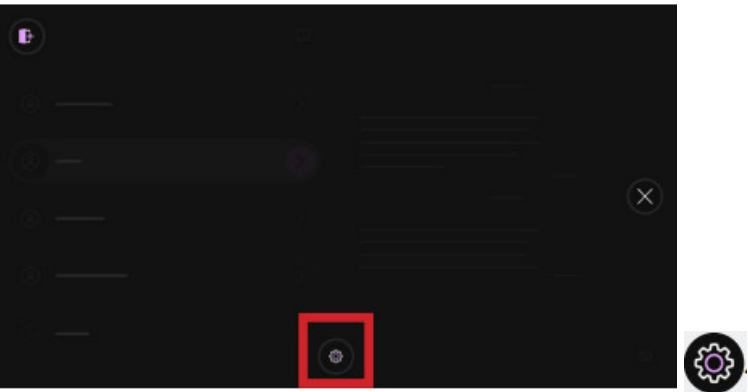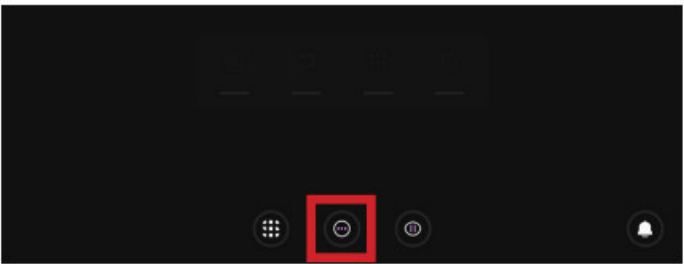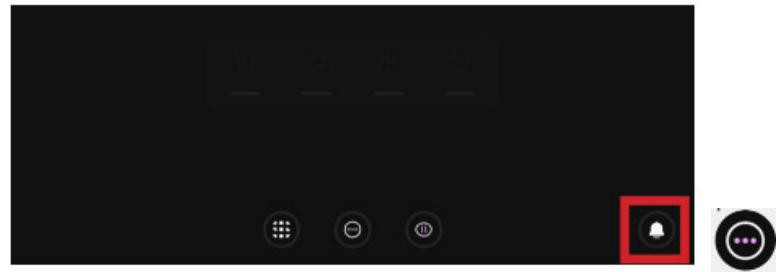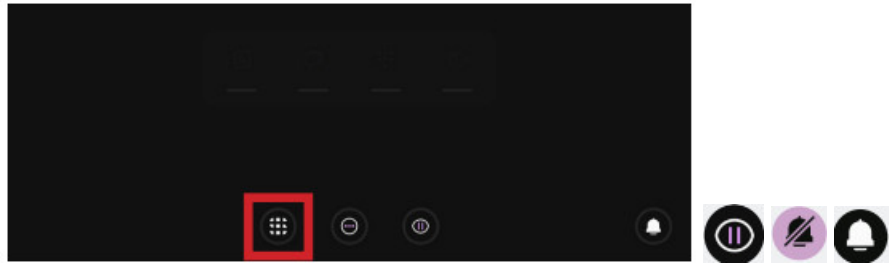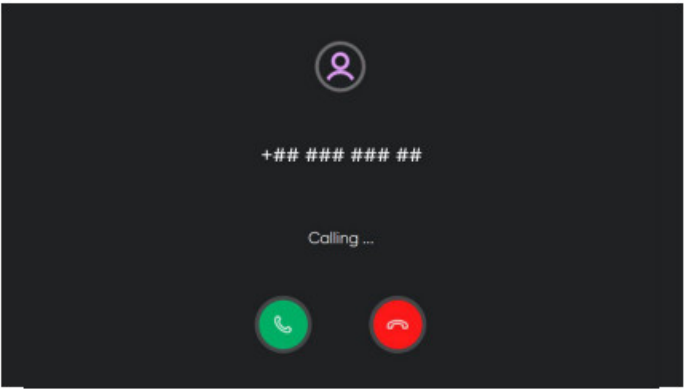

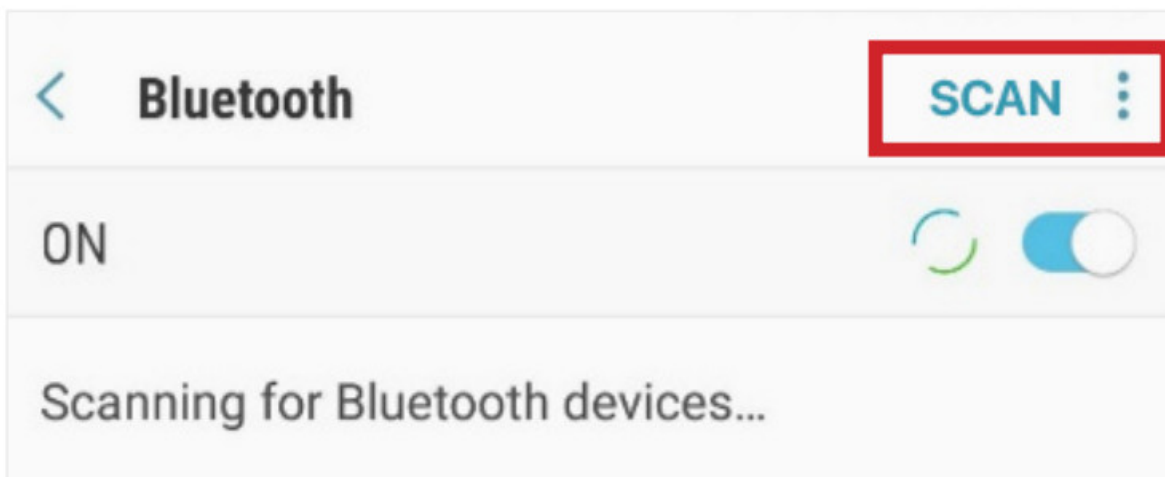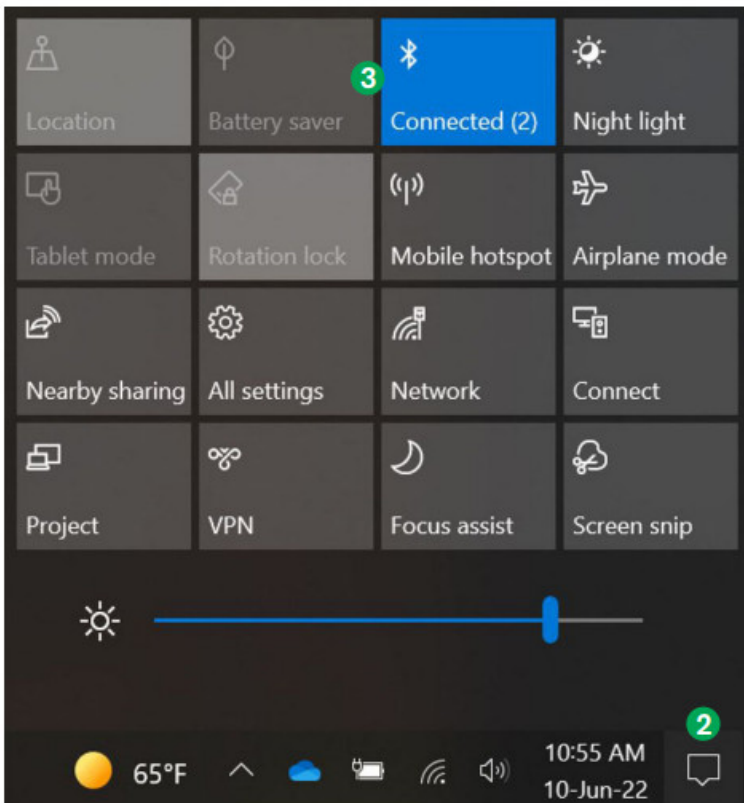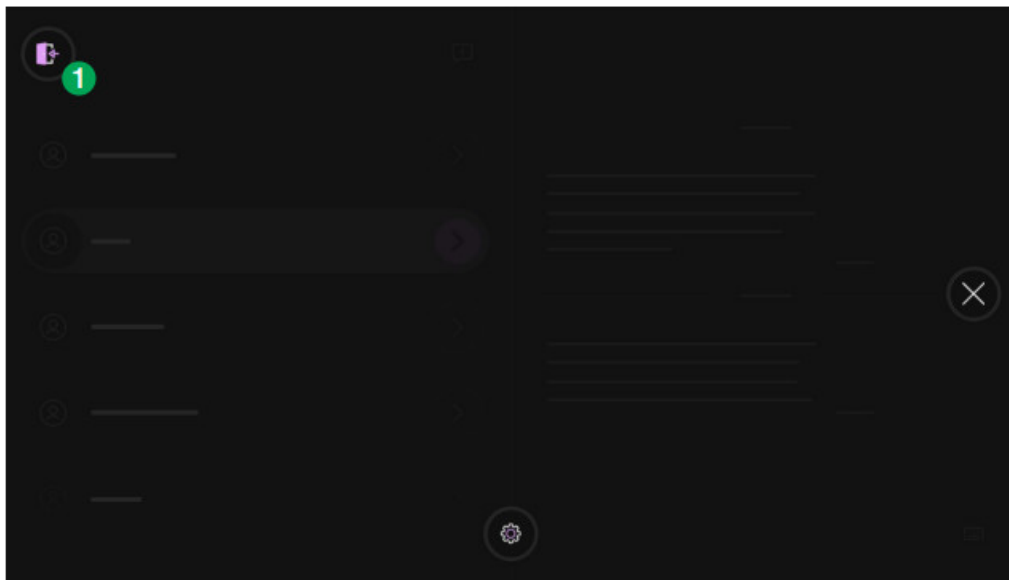




Documents / Resources

	tobii dynavox TD Phone Connection [pdf] Instructions I-13, I-16, TD Phone Connection, TD Phone, Connection
---	---

References

- [tobiidynavox](#)
- [Log into Facebook | Facebook](#)
- [User Manual](#)

[Manuals+](#), [Privacy Policy](#)

This website is an independent publication and is neither affiliated with nor endorsed by any of the trademark owners. The "Bluetooth®" word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. The "Wi-Fi®" word mark and logos are registered trademarks owned by the Wi-Fi Alliance. Any use of these marks on this website does not imply any affiliation with or endorsement.