




# Traceable 5665TR Three Channel Alarm Timer User Guide

[Home](#) » [TRACEABLE](#) » Traceable 5665TR Three Channel Alarm Timer User Guide 

## Traceable® 5665TR Three Channel Alarm Timer User Guide

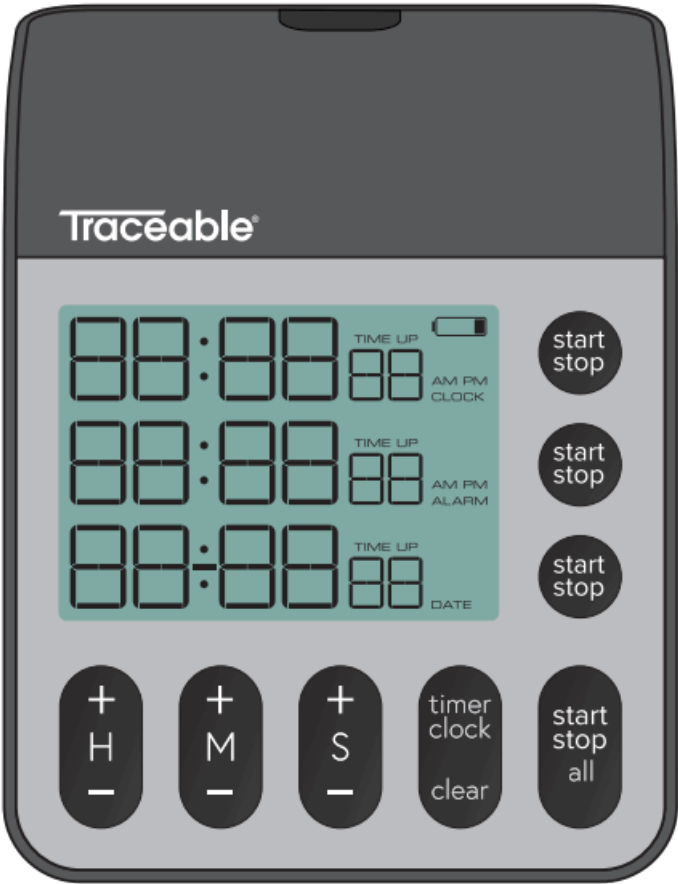


## Contents

- [1 5665TR Three Channel Alarm Timer](#)
- [2 INITIAL SETUP](#)
- [3 Documents / Resources](#)
  - [3.1 References](#)
- [4 Related Posts](#)


## 5665TR Three Channel Alarm Timer

Traceable® Products are ISO 9001:2018 Quality- Certified by DNV and ISO/IEC 17025:2017 accredited as a Calibration Laboratory by A2LA. Item No. 56000-13 Model No. 5665 ©2022 1065T8\_M\_92-5665-00 Rev. 0 09022022

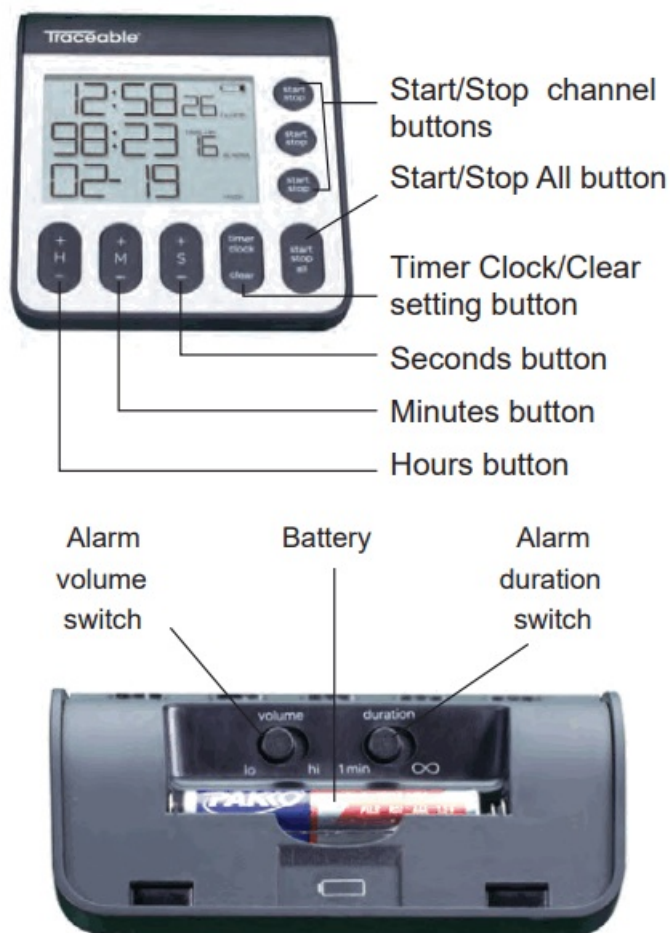


## INITIAL SETUP

1. Remove battery cover and install batteries.

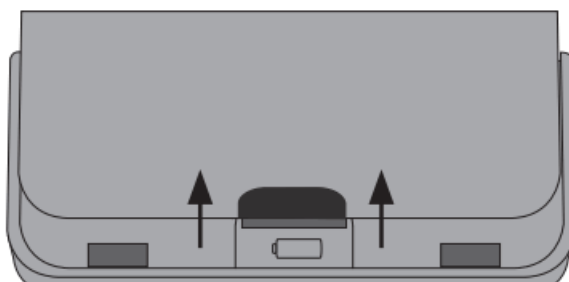

2. Use the volume slider to adjust buzzer volume.

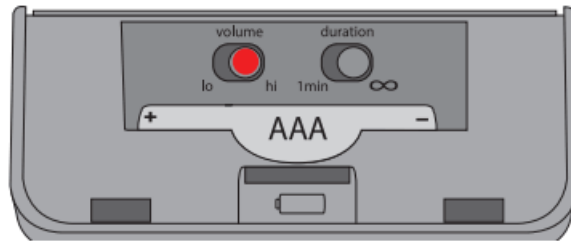


3. Use the duration slider to adjust alarm duration.




4. Replace battery cover.


5. Press TIMER/CLOCK until time appears on screen.

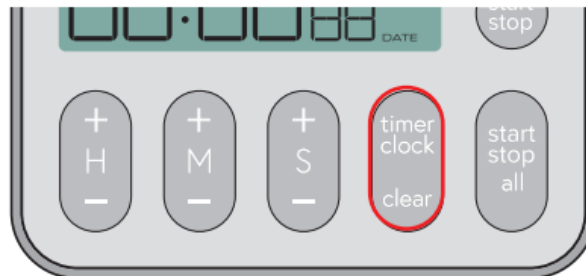


6. Hold TIMER/CLOCK until there is an audible beep.

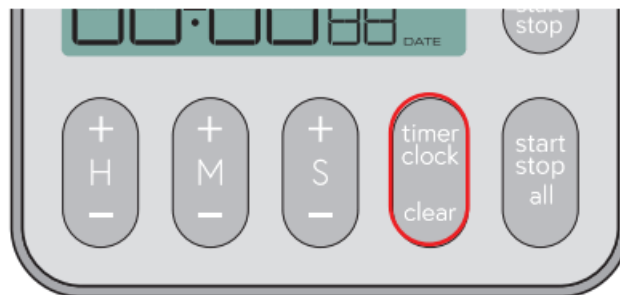

7. Use H/M/S buttons to set time; press TIMER/CLOCK when done.

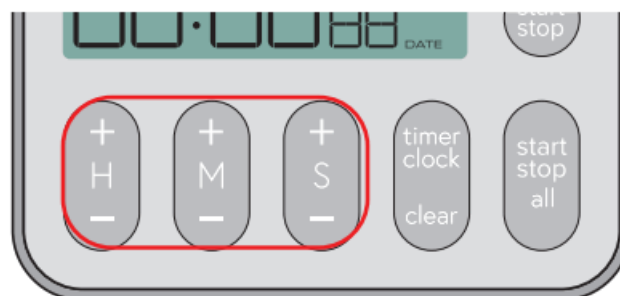

8. Use H/M/S buttons to configure alarm and date.

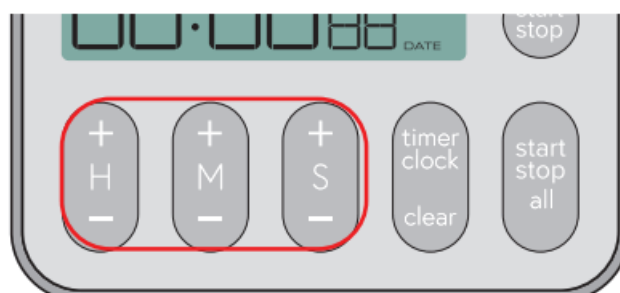

**Note:** To activate alarm, press START/ STOP ALL button when ALARM is lashing. ALARM will appear on display if alarm is active.


### SETTING A TIMER



1. Press TIMER/CLOCK until timer is displayed.
2. Press and hold START/STOP on the timing channel to be configured.
3. Use H/M/S buttons to configure timer.
4. Press START/STOP on timing channel to set timer.

**Note:** To recall memory, simply press START/STOP button after the timer has alarmed.

### STARTING A TIMER

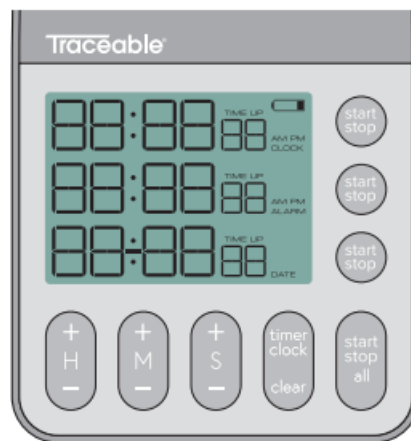

To start/stop an individual channel, press the START/STOP button for that channel. To start/stop all channels, press the START/STOP ALL button

**Note:** If a channel has not been configured with a set time, pressing START/ STOP on that channel or START/STOPALL will start a count up and TIME UP will flash on that channel.




# Traceable®

1.888.610.7664

[www.calcert.com](http://www.calcert.com)

[sales@calcert.com](mailto:sales@calcert.com)

**Documents / Resources**

  
  
  
Traceable® is a registered trademark of Traceable Systems, Inc. All other trademarks are the property of their respective owners.  
© 2015 Traceable Systems, Inc. All rights reserved.  
Model No. 5665TR  
Traceable Systems, Inc. 4000 West 10th Ave. Suite 100  
1-800-245-7664 • www.traceable.com • sales@traceable.com

## Traceable 5665TR Three Channel Alarm Timer [pdf] User Guide

5665TR Three Channel Alarm Timer, 5665TR, Three Channel Alarm Timer, Alarm Timer

References

-  [Calcert](#)

[Manuals+.](#)